



INFORMATIVE NOTE

Name: REFLEXX N80B

Date: 06/04/2021 Rev: 01

MANUFACTURER

REFLEXX SpA
Via Passeri, 2 46019 VIADANA (MN)
Tel.: + 39 0375 758891 Fax + 39 0375
464504 info@reflexx.com
www.reflexx.com



Pack		
Cod./Size	Box	Carton
N80B/S	8032891631065	8032891636060
N80B /M	8032891631072	8032891636077
N80B /L	8032891631089	8032891636084
N80B/XL	8032891631096	8032891636091

DESCRIPTION: ■ disposable glove powder free, not sterile ■ **Material:** nitrile ■ **Shape:** ambidextrous; five fingers, beaded cuff ■ **Colour:** blue ■ **Feature:** internally chlorinated; textured surface and textured fingertips ■ **Packing:** 100 pieces/box; 10 boxes/carton ■ **ISO 2859:** AQL 1.5 (level G1) freedom from holes; AQL 2.5 (level G1) major defects; AQL 4 (level G1) minor defects

CONFORMITY AND PERFORMANCES

Personal Protective Equipment cat III

0465 ANCI_CIMAC, Via Aguzzafame, 60/b 27029 Vigevano (PV)

Reg EU 2016/425 - EN ISO 374-1:2016 TYPE C - EN ISO 374-5:2016 - EN 374-2:2014 - EN 420:2003 + A1 2009

Medical Device

Reg EU 2017/745- EN 455 1,2,3 & 4

(CND T01020204 NON SURGICAL NITRILE GLOVES)



cat III
classe I



Reg CE 1935/2004 – Reg EU 10/2011 - DM 21/03/19. To such use carefully read the declaration of suitability for food contact (DICA).

DIMENSIONS		
Cod./Size	Length min(mm)	Width (±5 mm)
N80B/S (7)	240	85
N80B/M (8)	240	95
N80B/L (9)	240	106
N80B/XL (10)	240	116


Thickness mm	Palm (+/- 0.02)	0.06
--------------	-----------------	------

EN 455-2:2009/A1	Tensile strength	≥6 N
------------------	------------------	------

EN ISO 374-5:2016	Protection from microorganisms (bacteria and fungi)	Virus protection ISO 16604:2004
	pass	pass

EN ISO 374-1:2016				
Chemical compound :		EN 16523-1:2015 TYPE C		EN 374-4:2013
	COD	Performance level	Breakthrough time	% degradation
Sodium hydroxide 40%	K	6	>480'	6,5%

REG EU 2016/425

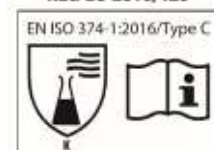

EN 420:2003 + A1:2009

RISKS FROM WHICH THE DPI IS DESIGNED TO PROTECT:

Chemical compound :	Protection
K Sodium hydroxide 40 %	CAS:1310-73-2 Protection from use of inorganic bases

RECOMENDATIONS OF USE Multi-purpose glove to protect the hand and/or the handled product, can also be used for general uses ■ The glove is suitable for contact with pure chemical products for a period of less than the permeation time indicated in the chemical resistance table ■ Results obtained under laboratory conditions may not reflect actual working time **WARNINGS** Before using the product, try it in the real condition of use ■ The glove is a device subject to wear, so frequent replacement is recommended and in all

cases where imperfections are observed ■ Wear gloves with dry, clean hands ■ In rare cases where transient hypersensitive skin reactions may occur, discontinue use and consult a doctor ■ The user must inform the manufacturer and the relative authority of any incidents occurring during the use of the protective equipment

STORAGE AND DISPOSAL Store the gloves in their original packaging in a dry place. ■ Do not expose to direct sunlight and heat sources. ■ The expiry date is 3 years from the date of manufacture indicated on the sales package, and is valid for the product properly stored ■ Gloves must be disposed of according to local regulations ■ Product Declaration of Conformity on www.reflexx.com